

Semi-automated coagulation analysers

▲BIO SOLE▲ 2 & ▲BIO SOLE▲ 4



www.biolabo.fr



BIO SOLEA 2 BIO SOLEA 4



Practical and compact

2 and 4 channel semi-automated analysers for haemostasis

Your coagulometer for routine Tests

- Opto-mechanical detection system
- Automatic start upon reagent dispensing
- Open system with no maintenance

2 channel coagulometer BIO SOLEA 2



4 channel coagulometer BIO SOLEA 4



BIOLABO
Diagnostics






All routine analysis

- PT (Prothrombin Time)
- APTT (Activated Partial Thromboplastin Time)
- TT (Thrombin Time)
- Fibrinogen
- Plus 4 free positions for coagulation factors



Simple and automatic

- Works with any standard pipette
- Self diagnostic when switched on
- TEMP LO message until the instrument reaches the working temperature
- Automatic incubation after insertion of sample cups
- Automatic start (measurement) upon reagent dispensing
- Results automatically calculated displayed and printed

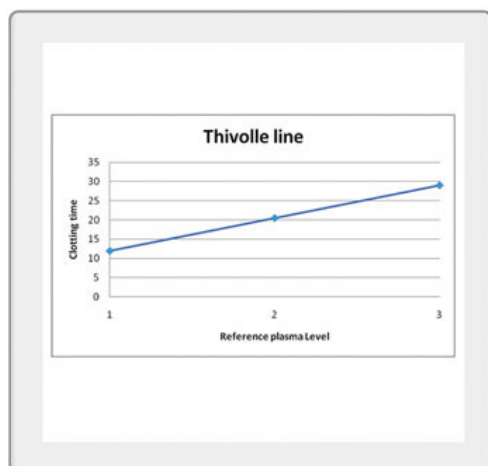

2 / 4 channel coagulometer BIO SOLEA 2 / 4 Technical features*

- Supplied with balls dispenser and a disposable starting pack
- Measuring method: opto-mechanical (optical and mechanical)
- Operating voltage: 110/220 V 50–60 HZ

BIO SOLEA 2

- Power consumption: 20 VA
- Net Dimensions: 31 cm x 12,5 cm x 11 cm
- Gross Dimensions: 40,5 cm x 30,5 cm x 18 cm
- Net Weight: 2,2 kg
- Gross Weight: 4,2 kg

BIO SOLEA 4

- Power consumption: 30 VA
- Net Dimensions: 33 cm x 25 cm x 11 cm
- Gross Dimensions: 40,5 cm x 31 cm x 27,5 cm
- Net Weight: 4,8 kg
- Gross Weight: 6,8 kg
- Magnetic stirrer for reagents

Printer (option)

Thermal printer
+ 1 Connecting CABLE + 2 rolls of thermopaper

*BIOLABO reserves the right to amend without notice the specifications of its products.



Commercial characteristics

- Thousands of units of these models had been sold world-wide under various and famous brand names.





Technical Service

- Clinical Chemistry specialists answer your technical questions. Please contact us for further details.
- Video tutorial available for an easy learning.



BIO SOLE 4

BIO SOLE 2



Les Hautes Rives
02160 MAIZY
France

Tél. + 33 (0)323 251 550
Fax + 33 (0)323 256 256

BIOLABO
Diagnostics




sales@biolabo.fr



support@biolabo.fr

Distributor:

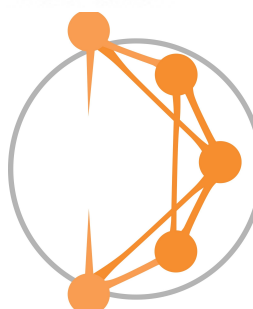

DIAGNOSTIKA PILIPINAS, INC.

5th Floor SBSI Bldg. 6023 Sacred Heart corner Kamagong Street, Brgy. San Antonio Makati City, 1203, Philippines

Main: +63 (2) 8638-0057
Logistics & Delivery: +63 (2) 8634-1487
Fax: +63 (2) 5637-7514
Service: +63 (9) 178161659

version: 22 / 01 / 2019

www.biolabo.fr